



VariFit™ – NIV mask



Critical Care ▪ Patient Interfaces



VariFit™ featuring AIR^ogel® technology

The VariFit is a single patient use hospital NIV mask featuring a cushion with AIR^ogel® technology for enhanced patient comfort. Colour coded elbows help the clinician to choose the correct mask for bilevel/CPAP using a single limb breathing system, or non-invasive ventilation using a dual limb breathing system. The breathable, quick release headgear is designed to fit easily to the patient and improve patient comfort.

Video available at www.intersurgical.nl



AIR^ogel®
Maximum comfort, minimum weight.
AIR^ogel® minimises pressure points and skin irritation

Advanced Cushion Technology™


The ergonomic cushion technology forms an effective seal for a comfortable fit

Active Headgear Connector™

Moves with the patient for a reliable seal

360° rotating elbow

For optimal tube positioning and easier mobility

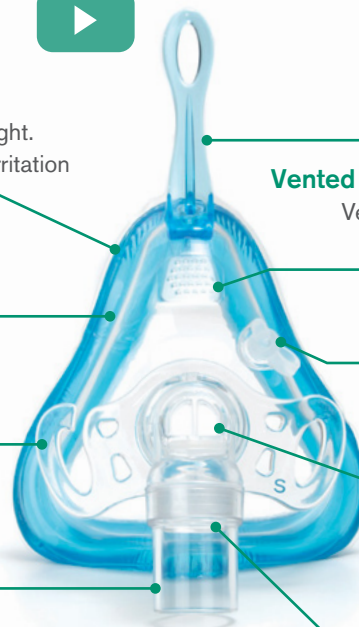

Flexible Spacebar™

The unique bendable design provides superior stability, maintaining a seal without contacting the skin

Vented Mask (2260000, 2261000, 2262000 and 2277000)

Vented mask option with an anti-asphyxiation valve for use with CPAP or bilevel devices

Oxygen Port

Allows easy delivery of supplementary oxygen

Anti-asphyxiation Valve

Prevents rebreathing in the case of inadequate ventilation

Colour Coded Elbows

To help the clinician choose the correct option

Code	Description	Colour coded elbow	Box Qty.
For use with single limb breathing systems that do not have a CO₂ leak port			
2260000	VariFit NIV vented mask with anti-asphyxiation valve, small adult	○	8
2261000	VariFit NIV vented mask with anti-asphyxiation valve, medium adult	○	8
2262000	VariFit NIV vented mask with anti-asphyxiation valve, large adult	○	8
2277000	VariFit NIV vented mask with anti-asphyxiation valve, extra large adult	○	8
For use with single limb breathing systems that have a CO₂ leak port			
2263000	VariFit NIV non-vented mask with anti-asphyxiation valve, small adult	◐	8
2264000	VariFit NIV non-vented mask with anti-asphyxiation valve, medium adult	◐	8
2269000	VariFit NIV non-vented mask with anti-asphyxiation valve, large adult	◐	8
2278000	VariFit NIV non-vented mask with anti-asphyxiation valve, extra large adult	◐	8
For use with dual limb breathing systems			
2265000	VariFit NIV non vented mask, small adult	●	8
2266000	VariFit NIV non vented mask, medium adult	●	8
2267000	VariFit NIV non vented mask, large adult	●	8
2279000	VariFit NIV non vented mask, extra large adult	●	8



Make an enquiry

AIRgel is a registered trademark of Sleepnet Corporation. Flexible Spacebar, Advanced Cushion Technology and Active Headgear Connector are trademarks of Sleepnet Corporation.

Please think before you print

Save energy and paper. If you must print this information sheet please print it double sided.

IS1.41 NL • Issue 6 02.19



Intersurgical Benelux B.V. Vluchttoord 28, 5406 XP, Uden, Nederland

T: +31 (0)413 243 860 F: +31 (0)413 243 869 info@intersurgical.nl www.intersurgical.nl

UK • Ireland • France • Germany • Spain • Portugal • Italy • Benelux • Sweden • Denmark • Lithuania • Russia • Czech Republic
Turkey • South Africa • China • Japan • Taiwan • Philippines • USA • Canada • Colombia • Australia