


STARMed CaStar R

Hood for non-invasive mechanical ventilation (NIV)

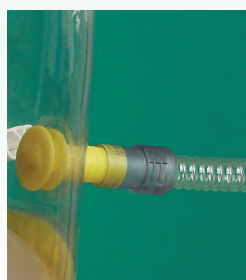

STARMed CaStar R – hood for NIV

StarMed CaStar R


The CaStar R is a comfortable, versatile and lightweight patient interface for non invasive mechanical ventilation (NIV) designed to offer an alternative to traditional methods of NIV delivery. The CaStar R is well tolerated for extended periods and can reduce the risks associated with both endotracheal intubation and face mask use. It can be used in both semi-recumbent and supine positions providing an ideal option for NIV delivery.

The CaStar R may be used for the treatment of a number of respiratory conditions including:



Pre and post extubation respiratory support, COPD exacerbation, thoracic trauma, asthma, post operative hypoxaemia, cardiogenic and non cardiogenic pulmonary oedema.

Technical Features

- 22M connectors compliant to ISO standard
- Internal volume 11 litres – weight 380g
- Comfortable and atraumatic underarm fastening straps
- Built in inflatable cuff
- Two sealed access ports for probes or catheters of 3.5-7.0mm Ø
- Patient access port with screw cap. Fitted with bidirectional anti-suffocation valve which opens automatically in the event of pressure failure
- Airtight sealing of the system is achieved by an elastic membrane making contact with the upper part of the torso, once the hood is inflated



Features and benefits

-   CE certified product
- Transparent**
For ease of patient monitoring
- Patient access port with anti-suffocation valve**
- Comfortable underarm fastening**

CaStar R hood with anti-suffocation valve and patient access port, ear plugs and measuring tape

Code	Size	Neck Circumference	Box qty.
CP211XXL/2RW	Extra Extra Large	≥50cm	5
CP211XL/2RW	Extra Large	45-52cm	5
CP211L/2RW	Large	40-47cm	5
CP211M/2RW	Medium	34-41cm	5
CP211S/2RW	Small	27-34cm	5
CP211XS/2RW	Extra Small	17-27cm	5

Single packed. Minimum purchasing package: 5 units. Custom-packs are supplied on specific request of the customer.

IS1.18_UK • Issue 3_05.13



Crane House, Molly Millars Lane,
Wokingham, Berkshire, RG41 2RZ
T. +44 (0)118 9656 300
F. +44 (0)118 9656 356

UK
info@intersurgical.co.uk
Italy
info@intersurgical.it

España
info@intersurgical-es.com
Portugal
info@intersurgical.pt

Lietuva
info@intersurgical.lt
Россия
info@intersurgical.ru

France
info@intersurgical.fr
South Africa
info@intersurgical.co.za

Japan
info@intersurgical.co.jp
Nederland
info@intersurgical.nl



Manufactured by: Intersurgical S.p.A.