



# Laryngeal mask airways

A comprehensive range of single use laryngeal mask airways



**Airway Management** • Airway Devices

## Laryngeal mask airways

A comprehensive range of single use laryngeal mask airways suitable for anaesthesia and emergency medicine, including the Solus™ Standard, the wire-reinforced Solus™ Flexible, the soft compound Solus™ Satin and Supra™ silicone options.

All Intersurgical laryngeal mask airways are single use, latex free and are supplied sterile, with the following features:

### Low friction material

Reduced insertion/extraction drag

### Classic cuff shape

Reliable consistent use

### Essential information

Visible on insertion for quick reference



### Integral inflation line\*

Reduces the risk of damage whilst in use

\* Integral inflation line refers to Solus products only


### High quality valve

To ensure continual cuff integrity

Video available at [www.intersurgical.com](http://www.intersurgical.com)

### Solus™ Standard

Using specially chosen materials, the Solus Standard is manufactured with the quality, flexibility and strength required to work effectively in a range of demanding clinical applications, including paediatric anaesthesia. Available in a full range, from adult to neonate.

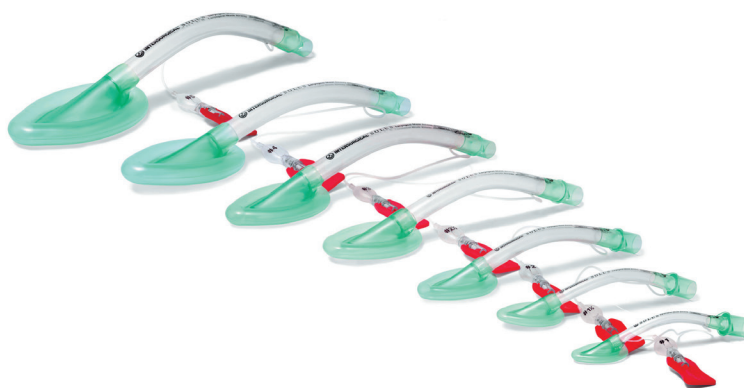

Code	Description	Size	Weight	Box Qty.
8005000	Solus Standard, large adult, laryngeal mask airway	5	70+kg	20
8004000	Solus Standard, medium adult, laryngeal mask airway	4	50–70kg	20
8003000	Solus Standard, small adult, laryngeal mask airway	3	30–50kg	20
8025000	Solus Standard, large paediatric, laryngeal mask airway	2.5	20–30kg	20
8002000	Solus Standard, small paediatric, laryngeal mask airway	2	10–20kg	20
8015000	Solus Standard, infant, laryngeal mask airway	1.5	5–10kg	20
8001000	Solus Standard, neonate, laryngeal mask airway	1	<5kg	20

All supplied sterile



## Solus™ Satin

Solus Satin is our range of single use laryngeal mask airways featuring a softer airway tube providing flexibility.



Code	Description	Size	Weight	Box Qty.
8005003	Solus Satin, large adult, laryngeal mask airway	5	70+kg	20
8004003	Solus Satin, medium adult, laryngeal mask airway	4	50–70kg	20
8003003	Solus Satin, small adult, laryngeal mask airway	3	30–50kg	20

## Solus™ Flexible

Solus Flexible provides an ideal solution for shared airway procedures such as ENT, dental, oro-maxillary and eye surgery where the rigidity of airway devices can obscure surgical access. The wire-reinforced tube permits flexion without kinking and can be moved at any time mid-procedure without concern for the cessation of gas flow.



Code	Description	Size	Weight	Box Qty.
8005001	Solus Flexible, large adult, wire-reinforced laryngeal mask airway	5	70+kg	20
8004001	Solus Flexible, medium adult, wire-reinforced laryngeal mask airway	4	50–70kg	20
8003001	Solus Flexible, small adult, wire-reinforced laryngeal mask airway	3	30–50kg	20
8025001	Solus Flexible, large paediatric, wire-reinforced laryngeal mask airway	2.5	20–30kg	20
8002001	Solus Flexible, small paediatric, wire-reinforced laryngeal mask airway	2	10–20kg	20

## Solus™ MRI-safe

A range of single use, sterile, laryngeal mask airways fitted with specially tested non-ferrous valves, guaranteed not to interfere with the magnet in an MRI scanner. The plastic valve has been selected to ensure full reliability throughout the shelf life of each device.



Code	Description	Size	Weight	Box Qty.
8005002	Solus MRI-safe, large adult, laryngeal mask airway	5	70+kg	20
8004002	Solus MRI-safe, medium adult, laryngeal mask airway	4	50–70kg	20
8003002	Solus MRI-safe, small adult, laryngeal mask airway	3	30–50kg	20
8025002	Solus MRI-safe, large paediatric, laryngeal mask airway	2.5	20–30kg	20
8002002	Solus MRI-safe, small paediatric, laryngeal mask airway	2	10–20kg	20
8015002	Solus MRI-safe, infant, laryngeal mask airway	1.5	5–10kg	20
8001002	Solus MRI-safe, neonate, laryngeal mask airway	1	<5kg	20

## Supra™ silicone

For those who require an alternative to reprocessable or PVC laryngeal mask airways, the Supra silicone range offers that. Available in seven sizes, from adult to neonate.



Code	Description	Size	Weight	Box Qty.
8105000 <i>NEW</i>	Supra, large adult, silicone laryngeal mask airway	5	70+kg	20
8104000 <i>NEW</i>	Supra, medium adult, silicone laryngeal mask airway	4	50–70kg	20
8103000 <i>NEW</i>	Supra, small adult, silicone laryngeal mask airway	3	30–50kg	20
8125000 <i>NEW</i>	Supra, large paediatric, silicone laryngeal mask airway	2.5	20–30kg	20
8102000 <i>NEW</i>	Supra, small paediatric, silicone laryngeal mask airway	2	10–20kg	20
8115000 <i>NEW</i>	Supra, infant, silicone, laryngeal mask airway	1.5	5–10kg	20
8101000 <i>NEW</i>	Supra, neonate, silicone laryngeal mask airway	1	<5kg	20

**S** All supplied sterile

IS6.1 NL • Issue 8 05.14



Intersurgical Ltd,  
Crane House, Molly Millars Lane,  
Wokingham, Berkshire, RG41 2RZ, UK  
T: +44 (0)118 9656 300 F: +44 (0)118 9656 356  
info@intersurgical.com www.intersurgical.com

**Deutschland**  
info@intersurgical.de

**France**  
info@intersurgical.fr

**Italia**  
info@intersurgical.it

**España**  
info@intersurgical-es.com

**Portugal**  
info@intersurgical.pt

**Nederland**  
info@intersurgical.nl

**Lietuva**  
info@intersurgical.lt

**Россия**  
info@intersurgical.ru

**Česká Republika**  
info@intersurgical.cz

**Taiwan**  
info@intersurgical.com.tw

**South Africa**  
info@intersurgical.co.za

**Philippines**  
info@intersurgical.ph

**Japan**  
info@intersurgical.co.jp

**USA**  
support@intersurgicalinc.com

**China**  
info@intersurgical-cn.com

